

## SE-M02 Beauty Easy With Multi-function



### TO BUYER:

Before linking or operating, please read this manual carefully. Thanks!

## User Manual

<http://www.SPAelf.com>

## **Summary**

The human body itself is the energy magnetic field, which is composed of many energy magnetic field, each part of the body needs have their own energy element, when the energy imbalance, can lead to the sub-health status. Traditional Chinese medicine says, fundamentality of su-health, that is short of yangqi (or energy), while many negative emotions (energy), in this highly competitive information age, the pressure of course, is inevitable. A comfortable stress-reducing method of their own is the physical and mental health.

Energy activation and conversion equipment has strong and stable energetic fields, easily blends with the human body magnetic energy, also can converse and enlarge the energy, transport the energy into the various organs, to purify the blood, lower blood viscosity, remove aging cutin of the heart, settle clutter thoughts, so that have good health and a good mood, free breathing fresh air, completely eliminate troubles which are brought from sub-health.

## **Working Theory**

Energy activation and conversion equipment integrates, bionic thermal, magnetic jade and ice technologies. Bionic thermal is conducted into the natural jade from Myanmar, rapidly destructs of disorder but beneficial trace elements in jade (zinc, magnesium, iron, copper, chromium, manganese, cobalt, etc.), makes them reorganize, to form 'Jade Energy Accumulator'. At the same time, bionic thermal effects on human cells and collagen, so skin cells quickly update to speed up the secretion of collagen by 300%. To form a special resonant magnetic field between 'Jade Energy Accumulator' and collagen, promote human body functions quicker in phase, to add necessary trace elements for the body, balance the imbalance of yin and yang, qi and blood, strengthen the skin self-absorption features.

## **Function**

**Care effect:** To promote infiltration and absorption of products, comprehensive care and treatment to enhance efficiency

**On Face:** skin scrubber, deep clean, penetration, restoration, anti aging.

**1, Skin Scrubber:** to removal the deadskin

**2, Deep clean:** to clean both superficial skin and deep skin completely

**3, Penetration:** Enhance skin absorption of products to make the skin red, white, glossy

**4, Restoration:** Thin pores, calm and repair telangiectasia, lessen thread veins.

**5, Antiaging:** Make up a large number of trace elements to skin collagen, promote the rearrangement of collagen material, in a short time to lift and tighten the cantus and drooping slack skin, lessen tiny wrinkle, dry wrinkle of the cantus. It takes effect 5 times better than normal.

**On Body:** Rapid charge to the cells tissue, a large number of supplementary energy for the body, the complete elimination of sub-health state

1, Purify the blood to remove waste, reduce blood viscosity, prevent high blood pressure and improve the anemia, and enhance immune function

2, To add trace elements necessary for the human body, add energy of jade 30 minutes, same as aerobic exercise outdoor 3 hours, to alleviate fatigue, the elimination of tension, the body more comfortable.

3, Menses adjustment and achesodyne, strengthen the liver and kidney function, and promote gastrointestinal peristalsis, regulate constipation

## **Market advantage**

1), To attract new customers and enhance the visibility of beauty salons

2), To enhance the beauty salon grade and expand profit space


3), To enhance the effectiveness of products, upgrade the product price, promote consumption.

4), To open new courses, new items.

5) , Can used together with any product, to enhance care efficiency, promote self-confidence of the beautician and then take the initiative to sell products.

6), One machine with multi function, very favorable, cater to customer's needs

## Control Panel

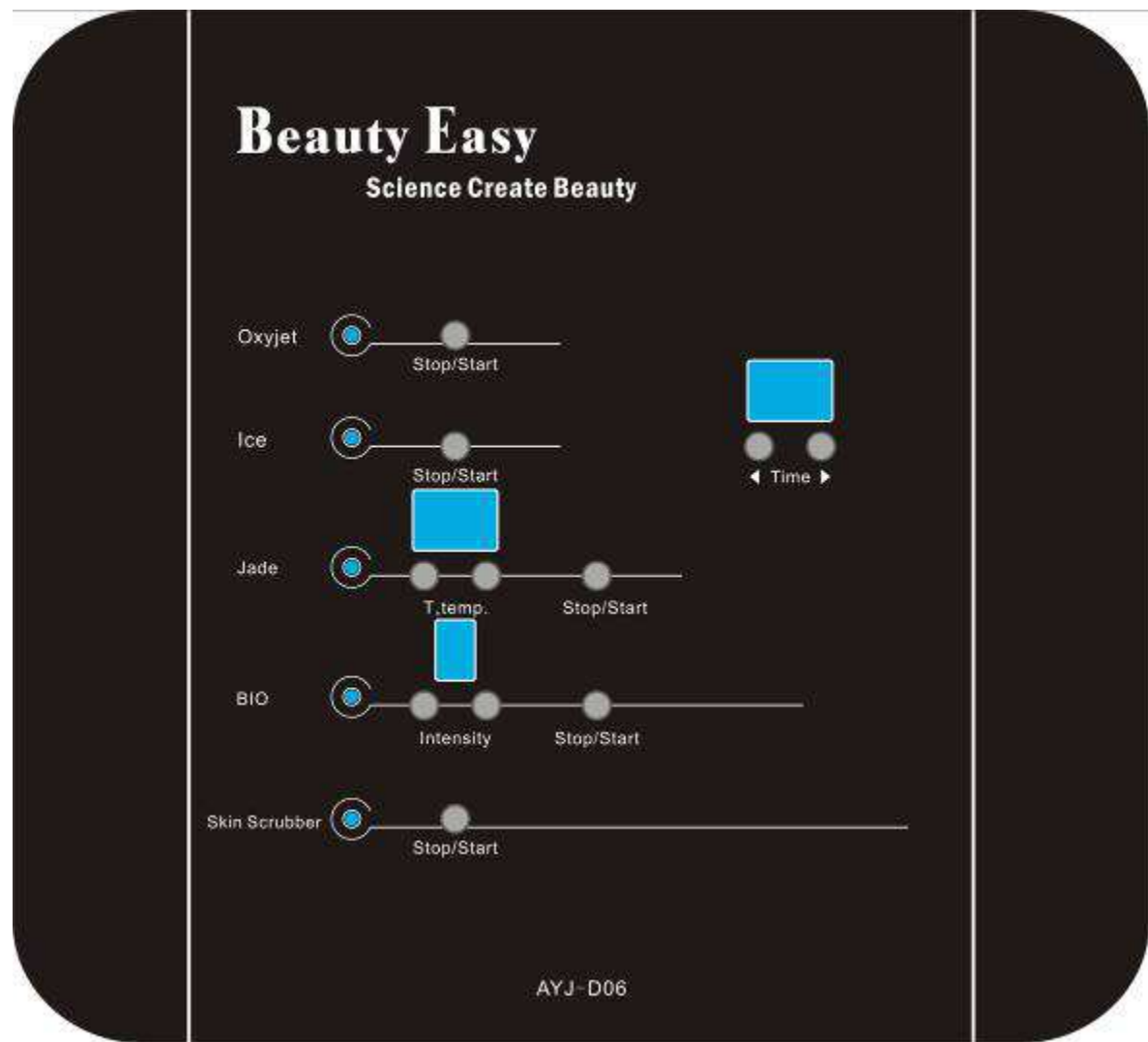

1. Oxyjet program: Stop/Start: To start or stop Oxyjet
2. Ice therapy Program: Stop/Start: To start or stop Cold hammer ()
3. Jade warm therapy Program: T. temp.: Turn down the temperature / Stop/Start: To start or stop Jade therapy.
4. BIO Magnetic energy Program: Stop/Start: To start or stop Magnetic therapy. Turn up the BIO intensity. Turn down the BIO magnetic current intensity.
7. Stop/Start: To start or stop skin scrubber
8. Display window for Time; Time (+); Time (-)

## Accessories Illustration



Strong activation rod



Jade head



3D ICE repair head



Skin scrubber



Oxyjet

### 1.Oxyjet

Inject the beauty product into the bottom of the epidermis via the micro-needle pressure oxygen injection. Which is to strengthen the cells, increase its resistibility, dispel the waste easily and promote its metabolism. So that it can treat and improve the problems such as dark and light skin, wrinkles, looseness, pigmentation and the conventional skin problems.

### 2.Strong activation rod+ Conductive probe:

To use bionic micro current, instantly open channels of pores, dredge the points, smoothen the meridian, rapidly pass life energy into deep cells tissue, to enhance cells permeability.

### 3. Jade head:

The pure natural jade is inlaid on the magnetic spin central, unique combination of 'magneto caloric effect' together with massage gesture, to carry through positive and negative energy conversion.

Within 3 seconds, re-organize the jade disorder trace elements, to form special 'Jade Energy Accumulator'. At the same time, bionic thermal effects on human cells and collagen, lets skin cells quick rejuvenate, speedily excrete 300% collagen.

### 4. 3D repair head

Through cold technologies, to calm and relax the skin, shrinken pores, penetrate nutrition into deep skin, lock inside water and nutrition, so as to get deep water supplement, water locking effect, promote absorption, oil controlling, skin tightening etc functions.

### 5. Skin Scrubber

Scrub the deadskin to make the skin fresh and the products can be easily to absorbed in inner skin.

## Accessories List

| No. | Description           | Qty. |
|-----|-----------------------|------|
| 1   | Strong activation rod | 2pc  |
| 2   | Jade head             | 1pc  |
| 3   | 3D Ice repair head    | 1pc  |
| 4   | Skin scrubber         | 1pc  |
| 5   | Oxyjet Gun            | 1pc  |
| 6   | Power cord            | 1pc  |
| 7   | English manual        | 1set |

## Specification

Voltage: 110V/60HZ,220V/50HZ

Power:150W

WoodenPacking:48X43X26cm

N.W.: 11KGS

G.W.:14.5KGS

# Operation Steps

**Applicable to the crowd:** all sub-health population

| <b>Special tip:</b> 1. Check all accessories whether they are connected right or not, clean or disinfect the accessories.<br>2. Plug in the power and switch on the general power on of the control panel.<br>3. Adjust the Time by Time button clicking according to treated course; |                                |                                                                                                                                                                                                                                                                                                                                                                                                                                                                           |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Step                                                                                                                                                                                                                                                                                  | Accessory                      | Operation course                                                                                                                                                                                                                                                                                                                                                                                                                                                          |
| 1 <sup>st</sup> step                                                                                                                                                                                                                                                                  | Skin scrubber                  | 1. Press stop/start button,<br>2. Time: 0-60 mins                                                                                                                                                                                                                                                                                                                                                                                                                         |
| 2 <sup>nd</sup> step                                                                                                                                                                                                                                                                  | Oxyjet                         | 1. Press stop/start button,<br>2. AC: 220V<br>3. 120W                                                                                                                                                                                                                                                                                                                                                                                                                     |
| 3 <sup>rd</sup> step                                                                                                                                                                                                                                                                  | Pairs of strong activation rod | 1. Press stop/start button under magnetic mode, then indicator flash, the client need to hold the conductive probe by hand.<br>2. Enter the course of 'opening points'. First press Intensity to adjust the current, then press stop/start, the right lower corner of the time, can see a small dot flash, that means the machine into working state.<br>3. Using the rod, along 10pcs of lymph drainage ways to do opening points, dredge meridians. <b>(10 minutes)</b> |
| 4 <sup>th</sup> step                                                                                                                                                                                                                                                                  | Jade head                      | 1. After finish, press stop/start to pause working, then press stop/start under Jade mode, the indicator flash.<br>2. Enter the course of energy conversion and deep purifying. According to the client skin quality, to press T,temp to adjust temperature, press stop/start to start working.<br>3. Use the jade head to transmit trace elements. <b>(20 minutes)</b>                                                                                                   |
| 5 <sup>th</sup> step                                                                                                                                                                                                                                                                  | 3D Ice repair head             | 1. After finish, press stop/start to pause working, under ICE mode, to press stop/start, then indicator flash, enter the course of 'Freezing therapy, resisting red veins floating'.<br>2. Press stop/start to start working, use 3d Ice repair head to calm and relax the skin. <b>(3 to 5 minutes)</b>                                                                                                                                                                  |

**Face:** Skin dry, short of water; wrinkles, sagging slack; red veins, large pores; spots pox;

**Body:** Early menopause, too much stress, insomnia, headache, memory loss, irritability upset, constipation, neck and shoulder pain, dizziness, gastrointestinal discomfort, neurasthenia, dysmenorrhea,

### **Conventional operation notes**

Before using the machine, please be sure to check whether the accessories are connected well or not.

1. First plug in the power, then switch on the power of the control panel, to avoid the phenomenon of reflux current.
2. Cables of accessories must be connected well with the sockets of the mainframe. Be sure to avoid inferior connection and lead to the useless of accessories
3. Be careful to move the machine.

Place it in the stable, dry place, when pushing it, please be sure to hold the handle, put the accessories onto the accessories holder, to avoid accessories falling to the ground and useless.

4. After finish, be sure to clean and disinfect the accessories. But don't use the alcohol or evaporating disinfectant to clean the mainframe. It will induce the mainframe broken.

Make sure to clean the residual products, to avoid products back to the mainframe, and damage the machine.

5. After starting up the machine, check if all intensity are set into zero.

During setting the intensity, please follow the customer skin quality, make sure it is comfortable. Do not set it too high to make skin with discomfort symptoms.

6. Strictly follow the time to operate the whole course.
7. Please keep treatment away of the skin with wounds. Do not operate it on Patients with facial paralysis.
8. Do not operate to whom with serious vasodilator.
9. Don not operate to whom with heart disease, high blood pressure, diabetic.
10. Before operation, please the customer and beautician take off all metal decorations.